



SUMMER IS COMING
The beach is waiting



**Skin is highly sensitive to
damaging environmental
factors like sun and sand**

A photograph of three young women standing on a sandy beach, seen from behind. They are wearing bikinis: black, white, and red. Their buttocks and thighs are covered in a thick layer of golden-brown powder or sand. The woman in the middle has her arms crossed over her hips. The woman on the right has her hands on her hips. The background shows the ocean and a clear blue sky.

**One of the most demanded
aesthetic treatments in summer
among young women is
BUTTOCKS TREATMENT**



Buttocks Treatment with JetPeel by TavTech


March 2024



Non-contact and non-invasive painless skin enhancing treatments which rejuvenates the skin and protects it from harmful environmental factors such as sand and sun



**The enduring
JetPro devices**

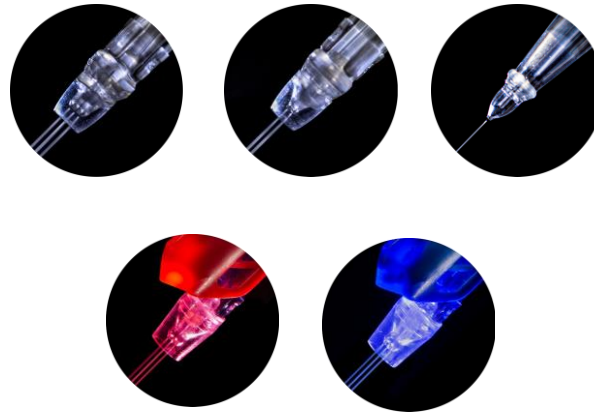


**The original patented
JetPeel handpieces**



**The premium
JetCare solutions**


Future Proof the Skin with JetCare Youth

The needle-free skin-enhancement JetPeel treatment is performed without touching the skin through a jet stream enriched with bio-active ingredients and a LED boost to stimulate the skin cells activities and to regenerate collagen

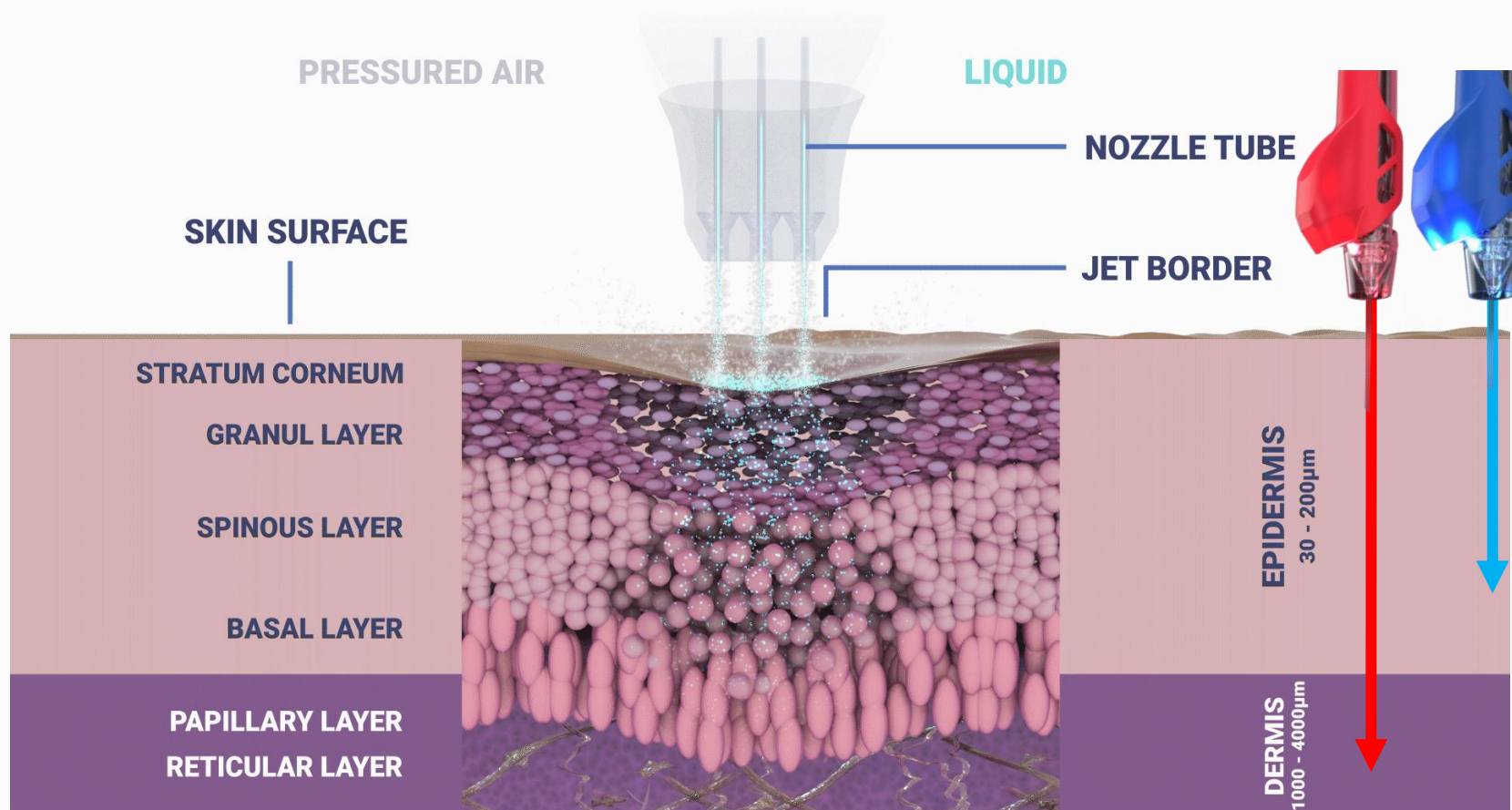

JetPeel by TavTech – Treatment Sequence

- Skin pressure massage to improve micro-circulation
- Pleasant exfoliation to remove dead skin cells for smoother and refreshed skin
- Trans-epidermal infusion of vitamins, nutrients and minerals
- Stimulate collagen and elastin production
- Slow down natural aging process
- Enhance skin elasticity
- Delay skin imperfections

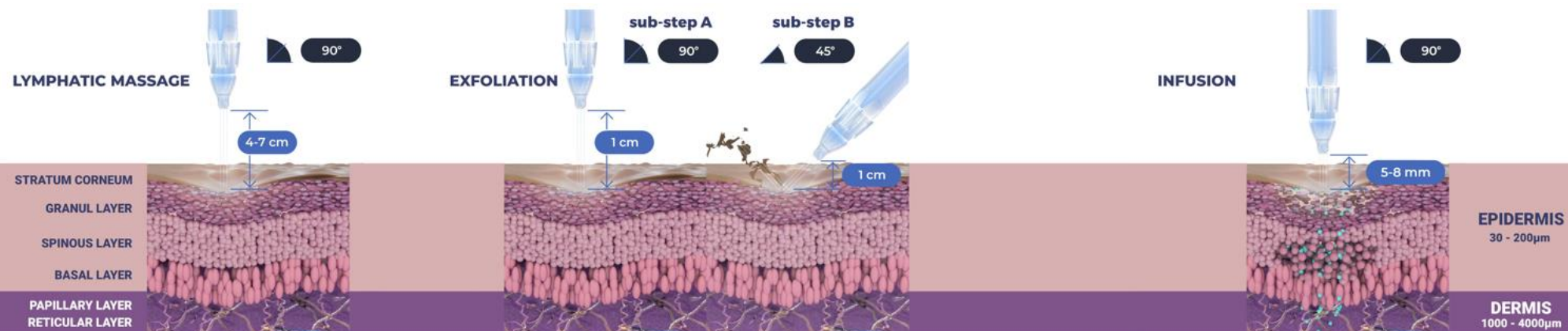

JetCare Youth enriched with Dead Sea minerals

is JetPeel's unique Skin Aging Prevention solution



Formulated with a blend of **21 Dead Sea minerals** and **natural antioxidants**, JetCare Youth comprises **3 different formulations:**

Dead Sea Hydro
Dead Sea Exfoliation
Dead Sea Nourishment



JetCare Youth modulates the biological function and prevents oxidative stress to maintain skin youth and vitality




This summer “PREJUVENATE” your buttocks with JetPeel by TavTech

- Stimulate collagen production and promote youthful appearance on the buttocks area
- Prevent loss of skin tone before the sun and sand accelerate the aging process
- Make the skin on the backside appear fuller, smoother and silky



Who is this treatment dedicated to?

Beach-loving bikini-wearers in their late teens, 20's and 30's to preserve the skin's natural youth

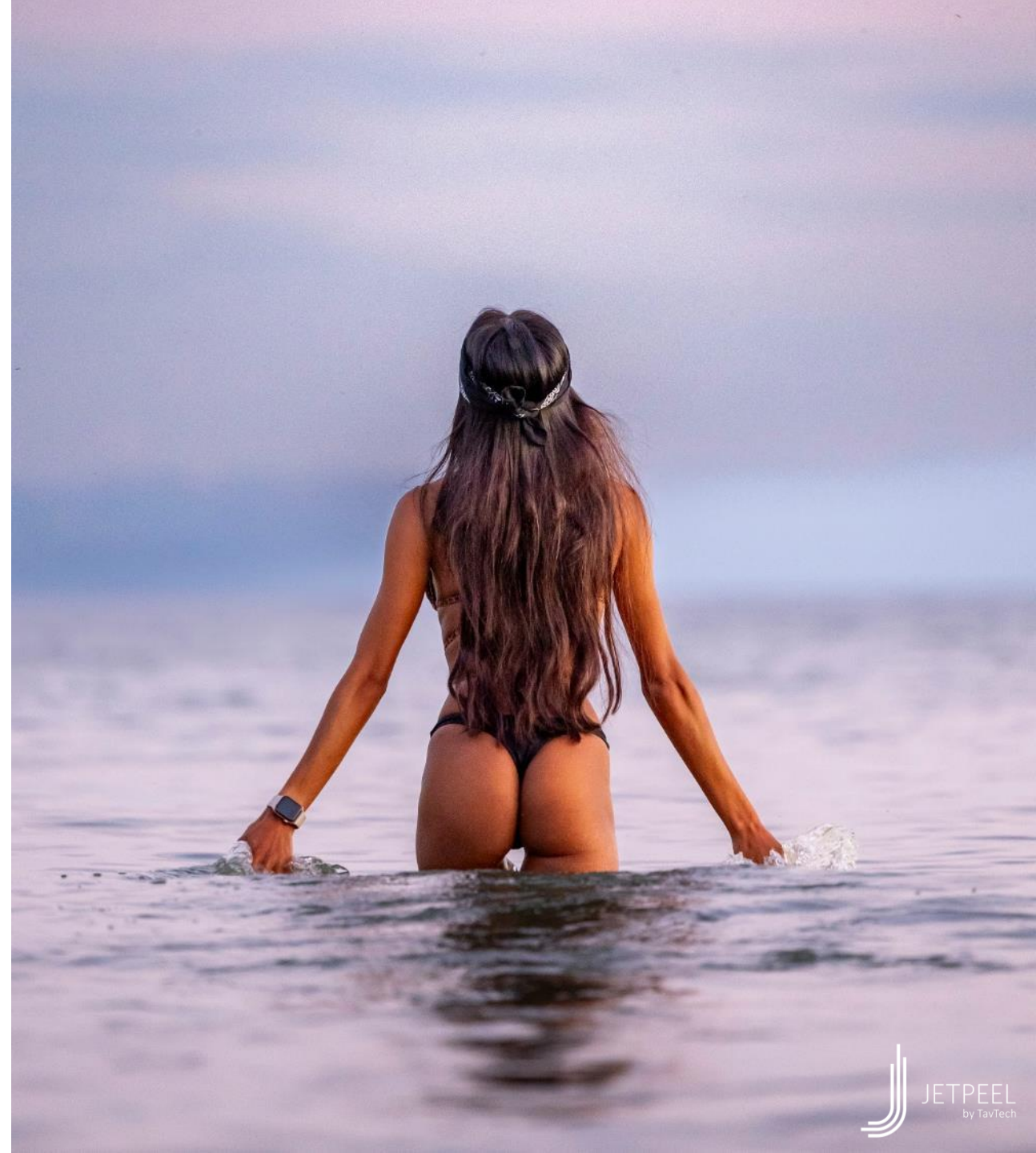

How JetPeel's prejuvenation buttocks treatment works?

Prejuvenate your butts in the bikini season **in just 3 steps!**

Treatment time - **30 minutes**



**STEP
1**

**LYMPHATIC
MASSAGE**

With Dead Sea Hydro



**STEP
2**

**GENTLE
EXFOLIATION**

With Dead Sea Exfoliation



**STEP
3**

INFUSION

With Dead Sea Nourishment



STEP 1- JetPeel 3D Massage

Stimulate micro-circulation and improve the skin tone

- **Enriched with Dead Sea mineral-rich salts:** Helps cleanse and detoxify the skin
- **Dunaliella HGL:** Helps reduce skin cell aging and protects from sun exposure damage
- **Organic Aloe Vera:** Antioxidant emollient that rapidly softens and calms the skin and provides antimicrobial action.



Use

10ml of the JetCare Dead Sea Hydro



Handpiece

TripleJet, held at 4 cm away from the skin at 90° to the skin



Time

8 minutes



STEP 2 - Exfoliation

Exfoliates, reactivates new skin cells, preserves a youthful skin appearance

- **Dunaliella HGL:** Helps reduce skin cell aging and protect from sun exposure damage
- **Allantoin:** Helps moisturize, exfoliate and smooth the skin surface
- **Lactic Acid:** Gentle AHA that helps remove dead cells, stimulates cell renewal, and boosts the skin's natural moisture
- **Organic Aloe Vera:** Antioxidant emollient that rapidly softens and calms the skin and provides antimicrobial action
- **D-Panthenol 75%:** Vitamin B₅, helps to seal in moisture, reinforce the skin's natural barrier, and protect from the effects of exposure and aging



Use



Handpiece



Time

Approximately 6 - 8 ml of Dead Sea Exfoliation in the sub-step A, leave the acid for 8 minutes and after that wash off with approximately 8ml of Dead Sea Hydro

TripleJet, held at 1 cm away from the skin and 90° to the skin in the sub-step A, and at 45° in the skin in sub-step B

14 minutes



Boost skin hydration and regeneration

- **Magnesium Sulfate:** Extra boost of this powerful anti-inflammatory mineral that helps remove excess sebum and toxins that contribute to skin damage.
- **Dunaliella HGL:** Helps reduce skin cell aging and protects from sun exposure damage.
- **Sepitonic™ M3:** Exclusive patented ingredient helps fight cell fatigue, boosts oxygen availability in the skin, and helps protect the dermis with anti-glycation action. It also acts as a carrier, helping to boost the effects of other anti-aging ingredients, strengthening the skin's resistance to damage from everyday elements.
- **Akoactive Queen:** Sourced from the fresh rose petal, this ingredient is a powerful skin hydrator, helping to reduce signs of skin fatigue, refreshing overall appearance and tone.



Use

Approximately 6 - 8ml of Dead Sea Nourishment



Handpiece

DoubleJet, held between 5 - 8mm away from the skin at 90° to the skin



Time

8 minutes

How many treatments should I have and how often?

This treatment, when performed with the **JetPro Duo** or with **JetPro ToGo**, simultaneously delivers **hydroporation** and **red LED photoporation therapy** to amplify the treatment results

6 treatment sessions to cover the bikini season

- **First 2 sessions:** twice a week, 3-4 days apart
- **3rd session:** after 8-12 days
- **4th session:** after 18-22 days
- **5th session:** after 28-32 days
- **6th session:** after 28-32 days


For best results, monthly treatments are highly recommended



- **Red LED light** affects mitochondrial function to produce ATP

**JetPro Duo and JetPro ToGo
provide the skin cells with:**

Nutrients (infusion) + light therapy
(photoporation) **for increased efficacy**

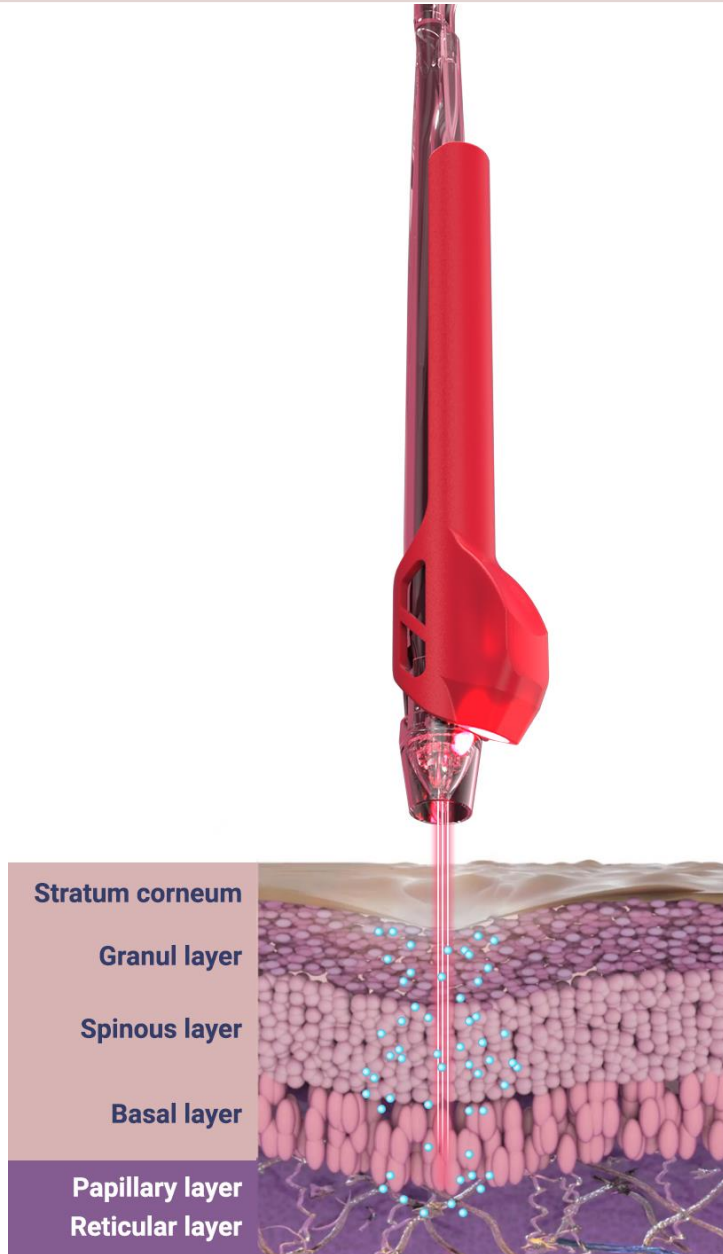



What results should we expect?

- Skin refinement and an even, radiant skin tone after just one treatment
- Skin protection from the damages of sun and sand exposure



**Jetpeel Buttocks Treatment
for healthy skin and self-confidence**



JETPEEL
by TavTech

Thank you!

